# **ARJOHUNTLEIGH**

**GETINGE GROUP** 

# MAXI SKY<sup>600</sup>

THE COMPLETE, CUSTOMISED CEILING LIFT SYSTEM



# A TOTAL CEILING LIFT CONCEPT

Maxi Sky<sup>600</sup> is a ceiling lift system that enables patients and residents to be transferred in a comfortable and dignified way. For the carer, it offers an easy-to-operate system that eliminates manual lifting.

A single carer can perform transfers quickly and efficiently, controlling all procedures with the handset. Maxi Sky<sup>600</sup> saves time, because it is always available for duty and within easy reach for the carer. And, when it's not in use, the unit takes up no precious floor space as it is discreetly out of sight at the recharging station on the overhead track.

What makes Maxi Sky<sup>600</sup> unique among ceiling lifts? It is a total concept. ArjoHuntleigh does much more than delivering an advanced ceiling lift, offering complete support every step of the way from track design to long-term system servicing.



#### · Designing the best system for your facility

We will help you create a layout that provides ideal coverage using our revolutionary modular track.

#### • Faster installation

Track installation is a quicker and easier process due to a new patented bracket design.

#### On track

The tight radius of the curved track sections enables layouts with superior space utilisation. Track sections fit together using a unique pin-locking system, and once in place, the seamless joints provide a smooth ride for the resident during lateral movement.

#### · Instant availability

The lift is always ready for use. After transfers, the lift can be programmed to return to the clip-on charging station, which can be positioned anywhere on the track.

#### Smooth operation and control

Regulation of lateral speed and spreader bar height is programmable from the handset. Electronic soft-start/ soft-stop enables the lift to be stopped smoothly and safely in the required position.

### • Simple maintenance

Battery changing is easy and requires only a screwdriver. A red light on the lift alerts the carer when servicing is required.

#### · The widest range of accessories

Access to our patient interface solutions - clip slings, loop slings, spreader bars and stretcher frame - gives you unrivalled flexibility to adapt to different needs.

### • Provision of staff training


We not only provide thorough guidance on how to use the system correctly, but also offer training on ergonomic working techniques.

### • Long-term system support

ArjoHuntleigh technical staff ensure that regular servicing and any unforeseen problems are dealt with quickly and efficiently.



# **ACCESSORIES** FOR OPTIMISED CARE

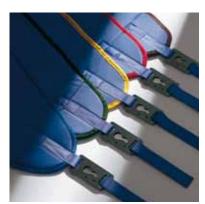

ArjoHuntleigh slings are colour-coded to clearly indicate sizes.



The 4-point Powered DPS spreader bar

Patient interface solutions are the key to optimising care for individuals. The comprehensive range of Maxi Sky600 accessories gives facilities flexibility to provide ideal solutions for different residents' requirements and to respond to changing needs. Sling solutions

A complete range of ArjoHuntleigh loop and clip slings covers the needs of all care facilities. Clip slings include standard padded and unpadded slings, special purpose slings for toileting, mesh slings for bathing and showering, and solutions for single and double amputees. The broad spectrum of sizes and sling types means that the ceiling lift can always be used with a sling that is optimised for the resident's comfort and safety.

### Spreader bars

ArjoHuntleigh loop and clip slings fit applicable 2-point spreader bars or 4-point DPS (Dynamic Positioning System) spreader bars. The DPS provides a safe and comfortable sling position for the patient during transfers and makes repositioning tasks easier for the carer. Powered DPS enhances DPS benefits by providing powered assistance to improve the precision and ease of resident repositioning. Using Powered DPS enables effortless, exact positioning of the resident at the press of a button.





#### **Walking Jacket**

Compatibility with the ArjoHuntleigh Walking Jacket means Maxi Sky<sup>600</sup> can also be used as an aid in rehabilitation. The Walking Jacket, which is applied to a 2-point loop spreader bar, provides safe and stable support during walking and stepping exercises.

#### Stretcher frame

Maxi Sky<sup>600</sup> is compatible with different stretcher solutions and is suitable for a number of types such as the scoop stretcher.

#### **Quick Connect**

At facilities where there is a regular need to switch between various spreader bars or to a stretcher, the Quick Connect accessory can provide a smooth, fast solution for changeovers.

#### **Electronic scale**

The optional Class 3-approved electronic scale allows residents to be weighed easily and efficiently during transfers.



The 4-point Flat DPS spreader bar, favourable e.g. low ceiling clearance.



Class 3-approved scale, used here on a 2-point spreader bar.





# A FLEXIBLE MODULAR TRACK SYSTEM



The exchanger and turntable enable a variety of different track solutions.



A turntable provides a choice of directions outside the resident's room.





continue on a different track after passing the exchanger.



Two, three or four tracks can be connected at various angles to the turntable, which is also powered by Maxi Sky<sup>600</sup>.



The gate, powered by Maxi Sky600, is a smooth way to combine a fixed track and an X-Y (traverse)

# **MAXI SKY IN CARE SITUATIONS**



If a resident has fallen, the long lifting strap (2.3 m) makes it easy and safe to perform a lift from the floor.



When full-room coverage is required, the X-Y (traverse) system is an ideal choice. With the unique traverse action, the track moves very easily with no risk of jamming. Here, it is used with an ArjoHuntleigh curtain to provide a high degree of privacy.




The track can be partially concealed in a suspended ceiling to provide a discreet solution.



There is a full range of Maxi Sky ceiling lifts that includes a portable lift, Maxi Sky<sup>440</sup>, with a safe working load of 200 kg (440 lbs) and, pictured here, Maxi Sky<sup>1000</sup>, which is designed to handle bariatric patients, and has a safe working load of 455 kg (1000 lbs).

## PRODUCT SPECIFICATIONS

#### Space requirement

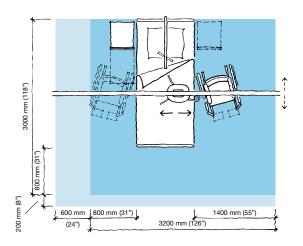

- Blue area shows the minimum working area required for the staff to be able to use the mechanical aids in an ergonomic way from one side.
- Light blue area shows required extension of working area to facilitate activities m either side to provide adequate access for the resident, mechanical aid and

The Maxi Sky<sup>600</sup> is part of the complete ArjoHuntleigh Maxi Sky ceiling lift range that includes the 1000 model, a bariatric ceiling lift, and the 440 model, a portable ceiling lift.

#### Mobility Gallery™

### Suitable for Doris and Emma



- D Doris, who has no capacity to support herself
- E Emma, who is almost completely bedridden and totally dependant Please contact ArjoHuntleigh for further information on the Mobility Gallery™

Product information	
Safe working load	272 kg (600 lbs.)
Unit weight (batteries included)	11.5 kg (28 lbs.)
Length	400 mm (15 ¾")
Width	256 mm (10 1/8")
Height	188 mm (7 ¾ ")

LED indicator when maintenance is required

Electronic microprocessor soft-start and stop motor control

Double centrifugal emergency brake system

Manual emergency lowering device (located on the motor cab)

Electrical up and down emergency buttons

Emergency stopping device (pull cord) accessible from the ground

Strap length up to 2.3 meter (90")

Lifting speed: 2.7 cm/sec. (1.2" / sec.) at 272 kg (600 lbs.)

3.5 cm/sec. (1.6" / sec.) at 200 kg (440 lbs.) 6 cm/sec. (2.3" / sec.) at 0 kg (0 lbs.)

Batteries: 2x5 Ah will average 150 cycles

(loaded at 75 kg / 165 lbs.)

Adjustable horizontal displacement speeds: 10, 15, 20, 25 cm/sec. (4, 6, 8, 10" / sec.). Preset at 20 cm/sec (8" / sec.)

User initiated return to charge station function with weight sensor cut-out of 9kg (20 lbs.) on/off programmable

ABS FR casing (fire retardant)

100-240 Vac / 50-60 Hz / 27 Va Charger unit

Clip on charger anywhere on the track

Class 2 insulated

Hand control Santopren tactile buttons / IP44

Infrared remote control (optional)

#### Approvals

CSA No 601.1, UL No 2601-1 certifications CE marked / ISO 10535

Only ArjoHuntleigh designed parts, which are designed specifically for the purpose, should be used on equipment and products supplied by ArjoHuntleigh As our policy is one of continuous development we reserve the right to modify designs and specifications without prior notice..

 ${}^{\otimes}$  and  ${}^{\top\!\!\!M}$  are trademarks of the ArjoHuntleigh group of companies  ${}^{\otimes}$  ArjoHuntleigh, 2010

# ARJOHUNTLEIGH **GETINGE GROUP**

#### www.ArjoHuntleigh.co.uk

#### ArjoHuntleigh UK - Head Office

310-312 Dallow Road, Luton.

Bedfordshire, LU1 1TD, United Kingdom Sales Phone: +44 (0) 1582 745700 +44 (0) 1582 745745 Sales Fax:

Rental 24hr Helpline Lo-call Phone: +44 (0) 8457 342000

+44 (0) 1582 745843 Rental Fax:

Email: sales.admin@huntleigh-healthcare.com GETINGE GROUP is a leading global provider of products and systems that contribute to quality enhancement and cost efficiency within healthcare and life sciences. We operate under the three brands of ArjoHuntleigh, GETINGE and MAQUET. ArjoHuntleigh focuses on patient mobility and wound management solutions. GETINGE provides solutions for infection control within healthcare and contamination prevention within life sciences. MAQUET specializes in solutions, therapies and products for surgical interventions and intensive care.

### ArjoHuntleigh UK - Service and Technical Centre

St Catherine Street

Gloucester, GL1 2SL, United Kingdom Phone: +44 (0) 8702 430430 +44 (0) 1452 525207 Service Fax: Email: ukservice@arjo.co.uk